


PROCEDURE FOR UPDATING THE USB CABLE DRIVER THROUGH THE INSTALLATION OF PROGUARD 4.7

Previous version of USB Cable driver: 4.40.1.0
Current version of the USB Cable driver: 6.0.0

Introduction

The new software PROGuard 4.7 has been developed to work with the most current version of the USB Cable driver (Version 6.0.0). This requires that after the installation of PROGuard 4.7, the new version of the USB Cable has to be installed as follows. Another option is to install the driver downloading it directly from the Contronics web site www.contronics.com



- 1 The PROGuard installation creates a directory called ProGuard drivers one level above Proguard as shown below. In this directory is the folder "USB Cable" that contains the driver installer with the name "CP210x_VCP_Win_XP_S2K3_Vista_7".


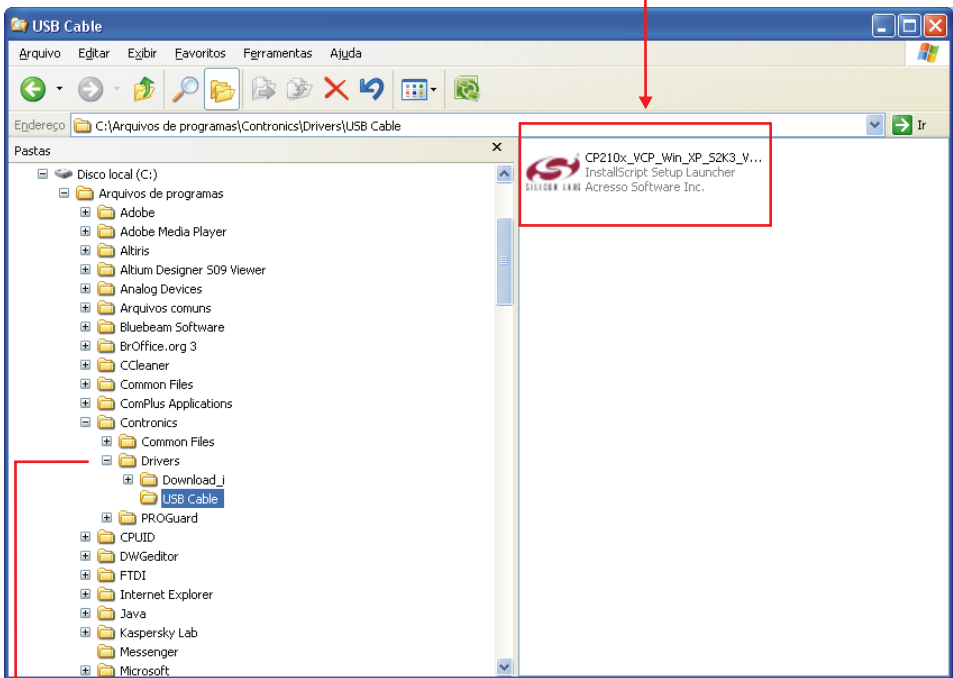



Directory created by PROGuard's installation, that contains the current drivers.

2 For the driver installation it is necessary that the device is disconnected from the computer. After unplugging it, run the installer "CP210x_VCP_Win_XP_S2K3_Vista_7"

3 Follow the steps:

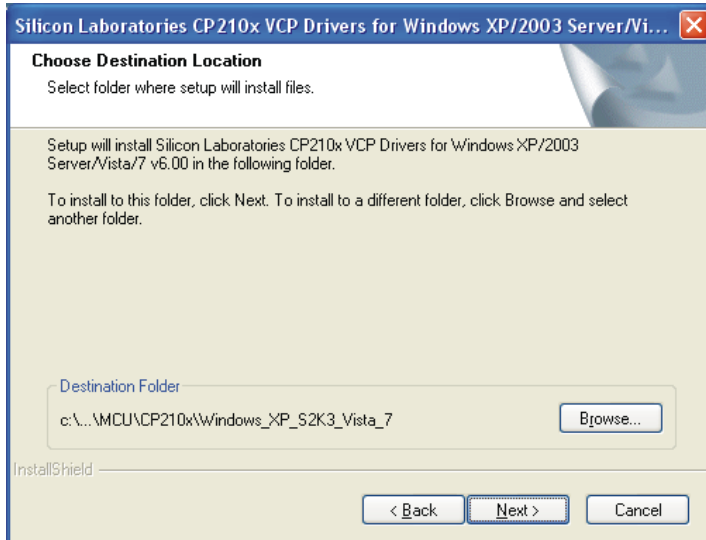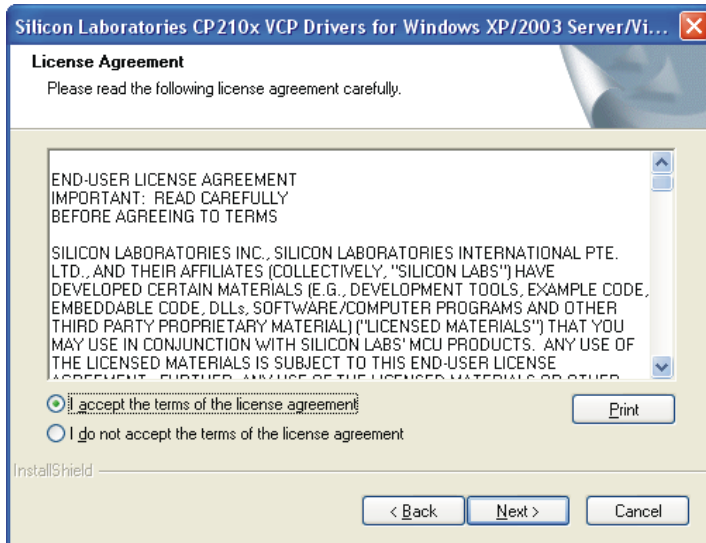

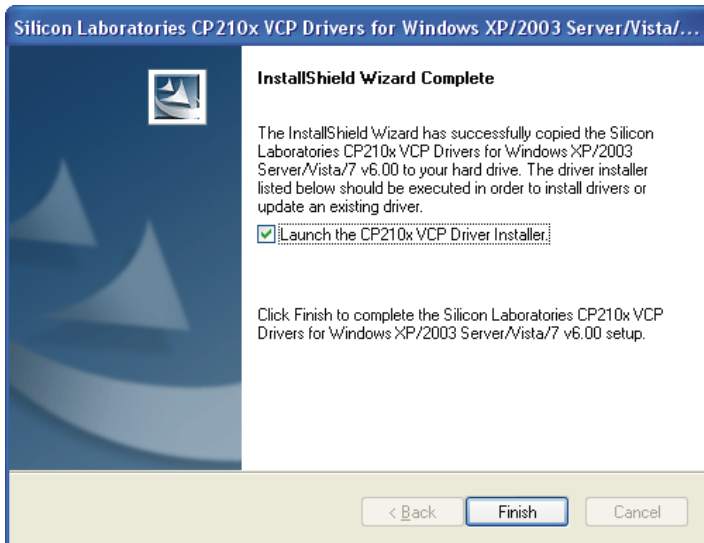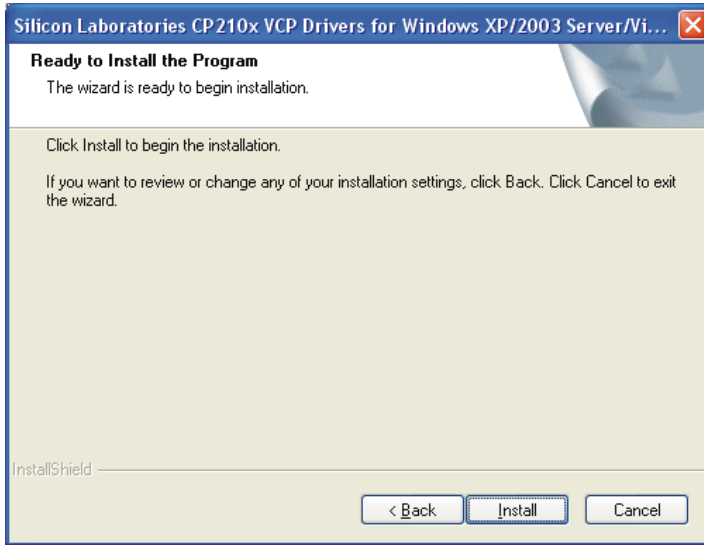

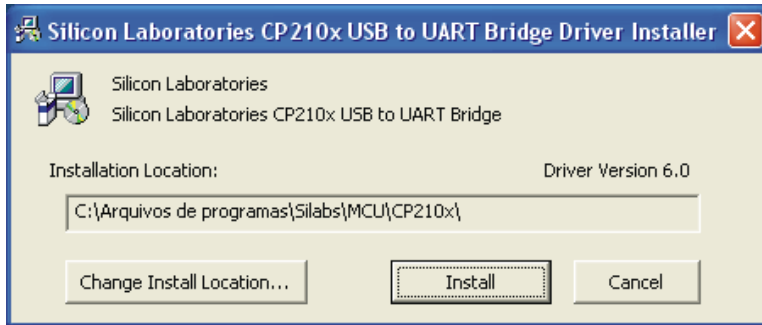

The installation is finished. If necessary restart your computer for installation to take effect.



www.contronics.com

Contronics Technologies INC.

8121 NW 68th Street
Miami, FL 33166 · USA
Phone: (305) 994-3535